

PF SERIES

Exam Table

- Choice of three different tops
- Beautiful wood base with choice of finishes
- Up to 4 powered motions



PF400
shown in Seafoam fabric
with optional Ebony
Lacquer and Locking
Casters



PF200



PF250




PF400



ISO 13485:2003
FM 584901

made in the USA
with US & imported parts

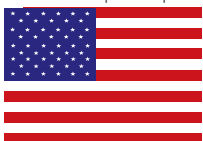

OAKWORKS[®]
MEDICAL
Position for Success

PF SERIES

Exam Table



ZERO|GRAVITY

Electric 4 Section Top



MODEL PF400

#70882-T* 31" W 4" Comfort
#73309-T* 29" W 4" Comfort

Motion 1	Electric Height: 22"-38" (56-97cm)
Motion 2	Electric Tilt: 20°/20°
Motion 3	Electric Backrest: 0-53°
Motion 4	Electric Leg rest



ZERO|GRAVITY

Lift-Assist Backrest Top



MODEL PF250

#70889-T* 31" W 3" Comfort
#73409-T* 31" W 4" Comfort
#73454-T* 29" W 3" Comfort
#73566-T* 29" W 4" Comfort
#73583-T* 27" W 3" Comfort
#73584-T* 27" W 4" Comfort

Motion 1	Electric Height: 22"-38" (56-97cm)
Motion 2	Electric Tilt: 20°/20°
Motion 3	Lift Assist Backrest: 0-80°

Flat Top



MODEL PF200

#70890-T* 31" W 3" Comfort
#73410-T* 31" W 4" Comfort
#73567-T* 29" W 3" Comfort
#73568-T* 29" W 4" Comfort
#73585-T* 27" W 3" Comfort
#73586-T* 27" W 4" Comfort

Motion 1	Electric Height: 21"-37" (53-94cm)
Motion 2	Electric Tilt: 20°/20°

PF SERIES EXAM TABLES

FEATURES	
Top Dimensions	see above for widths x 73" (185cm) long
Foot Control	Controls all powered movements
Weight Capacity	550 lb (250kg)
Product Weight	235 lb (106kg)
Padding - PF400	4" (10cm) Comfort Foam™
Padding - PF250 & PF200	3" (8cm) or 4" (10cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Warranty	3 years parts, 2 year labor
Safety Listings	FDA listed, CE marked, cETLus approved, CB Certificate

OPTIONS <i>must be ordered with table.</i>		#
4" (10cm) Locking Casters	To easily move your table in the room	69926
Wood Finishes	Choice of 4 stains or 4 lacquer finishes	73136
ACCESSORIES <i>order any time.</i>		#
Arm Rests (PF400)	Auto-adjusting & removable for easy access	74734-T*
Universal Table Extender with 4" Padding	21" (53cm) W x 12.5" (32cm) L	74736-T*
QuickLock™ Head Rest (black) & AeroCel Pad	Comfortable, adjustable prone or supine support.	PKG6914-T*
Paper Roll Holder	Vanilla - Accommodates 21" (53cm) roll (roll not included)	11580
Paper Roll Holder	Black - Accommodates 21" (53cm) roll (roll not included)	11349

*Wood finish and fabric color code must be specified.

The PF Series Exam Tables are designed to maximize comfort for extended treatments while offering the beauty of wood to create a Spa-like atmosphere. The fully electric design comes with a super soft and supportive padding system.

CASTERS #69926



OPTIONAL

HAND CONTROL

Add a hand control to control all functions.



enquire for
part number

OPTIONAL

ADJUSTABLE SIDE ARM RESTS #66760-T*

These incredibly versatile, removable arm rests rotate horizontally and always stay level. 7" (18cm) wide at widest point x 23.5" (60cm) long



(PF 400 Only)

ACCESSORY

PAPER ROLL HOLDER

#11580 vanilla #11349 black

This simple design is adjustable to accommodate different size rolls up to 21" (53cm) wide & plugs into the table's outlet holes. *(paper roll not included)*



ACCESSORY

QUICKLOCK™ HEAD REST & AERO·CEL PAD

#PKG4186-T*


The platform allows precise cervical flexion and height positioning over a wide range for both prone and supine treatments.



ACCESSORY


UNIVERSAL TABLE EXTENDER #68697-T* 3" padding #68735-T* 4" padding

Perfect as a head rest or foot rest, it has two rods that fit into outlets on the table. Extends table to 84". 21" (53cm) W x 12.5" (32cm) L

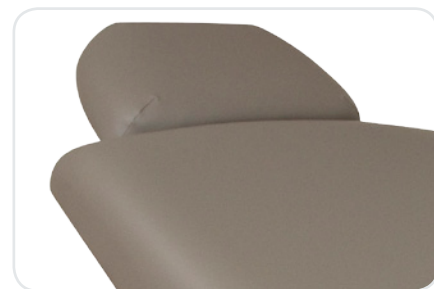


ACCESSORY

PADDING

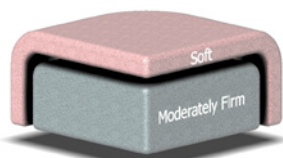

3" Comfort Foam™

More comfort and support
(Not available on PF 400)

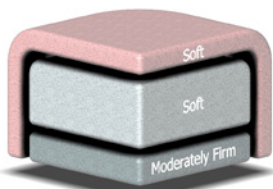

4" Comfort Foam™


Superior comfort and support

STAINS & LACQUERS


Oakworks[®] finishes are highly water resistant and can be wiped down with disinfectants as needed. Stains allow the natural grain of the wood to show through and create depth. Lacquers are more opaque and have a more solid appearance.

Standard
Natural
Finish

Optional Stains

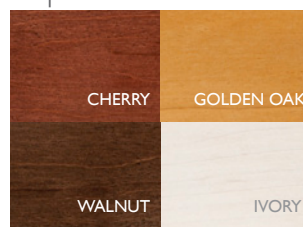


Optional Lacquers

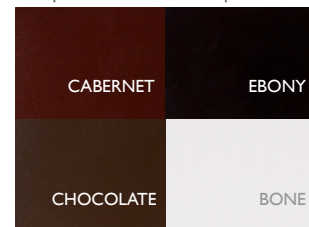

TERRATOUGH™ FABRIC

Meets the requirements of ISO 10993-1:2009 Biological Evaluation of Medical Tables and passes testing for Cytotoxicity, Skin Irritation and Sensitization as well as CA TB117 Fire Retardance testing. For disinfecting purposes you may use Protex, MadaCide, MadaWipe, Accell TB, Virox[®] (Canada), PDI Duper Sani-Cloth Germicidal wipe or a diluted 1/10 bleach solution.

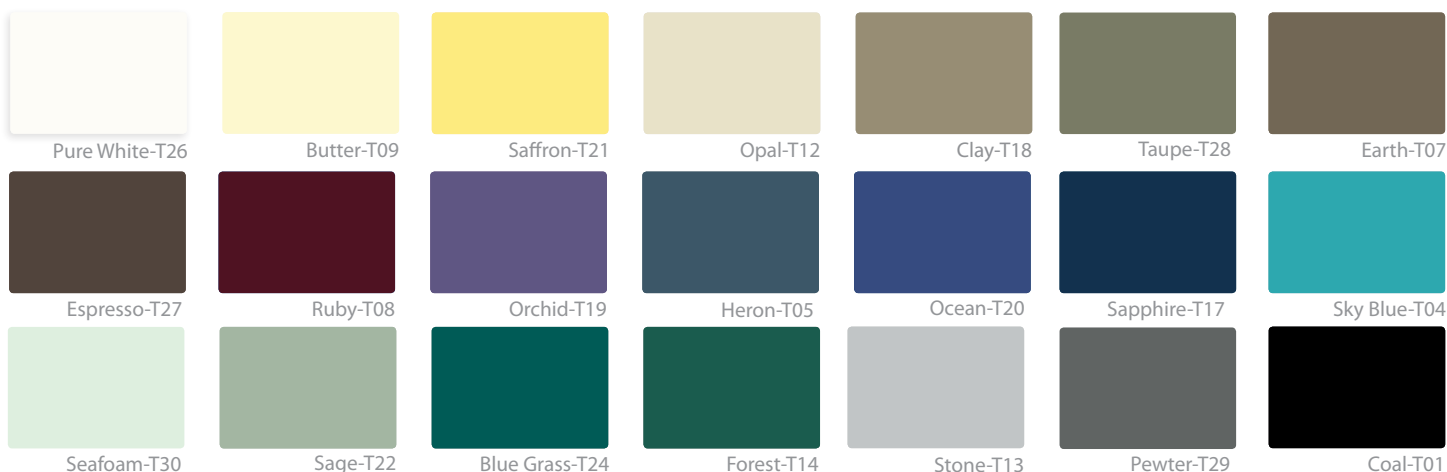

We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.

CONTACT:

*View our full line of
Medical products*



Scan this code with
your smart phone

#72703

